


ADVANCIS

Silflex Silicone Contact Layer Dressings

Silflex is an atraumatic wound contact layer dressing made from polyester mesh coated on both sides with Silfix soft silicone.

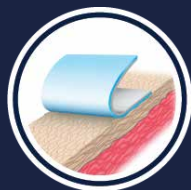

Atraumatic Silicone Coating

Reduced pain during application and dressing changes



Non-adherent

Does not adhere to the wound bed or secondary dressing



High Exudate Transfer

Large open pores allow exudate to easily pass through to the secondary dressing



Easy to Use

Easy to handle and does not adhere to itself



Conformable

Highly conformable for increased patient compliance¹²



NPWT Compatible

Can be used in conjunction with Negative Pressure Wound Therapy¹³



ADVANCIS

SILFLEX

Soft silicone wound contact layer



Description

Silflex is a wound contact, non-adherent dressing made from polyester mesh coated with an atraumatic, hydrophobic soft silicone, making it useful for patients with very thin, friable tissue.

Silflex is designed to gently adhere to the skin surrounding a wound and not to

the wound bed to reduce the pain and trauma associated with dressing change. With 2mm pore holes Silflex allows the passage of exudate through the dressing.

Due to its gentle application and removal, Silflex is an ideal solution for treating and preventing skin tears.

Ordering Information

Stock Codes	Size	Pack Size	NHSSC Code	PIP Code
CR3922	5cm x 7cm	10	EKH028	338-4385
CR3923	8cm x 10cm	10	EKH029	338-4393
CR3924	12cm x 15cm	10	EKH030	338-4401
CR3925	20cm x 30cm	10	EKH032	338-4419
CR4006	35cm x 60cm	5	EKH096	348-2205

Suitable for Wound Types:



ADVANCIS

SILTAPE

Soft silicone fixation tape



Description

Siltape is a soft silicone fixation tape delivered on a roll with perforations for ease of use.

Siltape is made from Silicone that is gentler on removal unlike other adhesives

that may cause tissue damage when removed and is particularly useful in patients with very thin, friable skin.

Indications

Siltape can be used in any situation where a regular adhesive tape would be used where friable or sensitive skin is an issue.

- For taping down eyelids for theatre.
- For securing fistula needles as it can be repositioned.
- For protecting friable tissue from pressure damage or sores from tubing; this is especially useful on the face and tops of ears. Siltape can be applied directly to the skin and tubing placed on top.
- For scar management

Ordering Information

Stock Codes	Size	Pack Size	NHSSC Code	PIP Code
CR3938	2cm x 3m	12	EHR053	344-0336
CR3939	4cm x 1.5m	12	EHR054	344-0328

Suitable for Wound Types:



ADVANCIS

WOUND PROBE

15cm sterile wound size assessment probe



Description

The Advancis Wound Probe is a single-use, disposable device designed for wound assessment.

It is used to measure the dimensions of a wound, including its width, length, and depth, to provide a comprehensive evaluation of wound size.

The probe body is strong, long, narrow and flexible, allowing for the assessment of single tracked sinus wounds, where only a primary track is present.

The probe can also be used to gently pack dressings.

To reduce trauma during assessment, the probe is smooth and has no raised or indented markings. The ergonomic probe handle is ridged and designed to be easily gripped between finger and thumb while wearing medical gloves.

Ordering Information

Stock Codes	Size	Pack Size	NHSSC Code
CR4223	15cm	10	ELZ864
CR4236	15cm	50	ELZ865

Please note - this product is not available for sale in Australia.

Suitable for Wound Types:



ADVANCIS

ADVASIL CONFORM

Silicone gel sheet for scar management



Description

Advasil Conform is a sheet of medical grade, self-adherent silicone gel for the management of scar tissue or recently healed wounds. Silicone gel sheets are reported to help reduce the effects of scarring if worn for prolonged periods. Hypertrophic and keloid scarring may particularly benefit from Advasil Conform.

Composed of a conformable layer of translucent Silfix silicone gel approved for skin contact, and a smooth polyurethane outer layer, Advasil Conform reduces pain and trauma at dressing changes, improving wear ability and patient compliance.

Advasil Conform is supplied sterile for improved hygiene, it is also washable and can be used for up to 28 days once open.

Indications

For the management of hypertrophic scars and keloid scars.

Ordering Information

Stock Codes	Size	Pack Size	NHSSC Code	PIP Code
CR3845	10cm x 10cm	5	EKH006	301-1640
CR3832	10cm x 15cm	5	ELY178	329-5979

Suitable for Wound Types:

